



KONICA MINOLTA

ARE You Satisfied with Your Current Print Quality & Cost




Introducing DNP Range FOTOLUSIO Printers

The first thing to consider when thinking about FOTLUSIO prints for your business is its outstanding print quality which helps you offer prints to match your reputation. Then consider FOTOLUSIO printers for its seamless Network print capabilities to connect to all your imaging modalities not to target that FOTOLUSIO printers offer the highest print quality at an affordable cost.

Ultrasound Imaging




PET CT Imaging

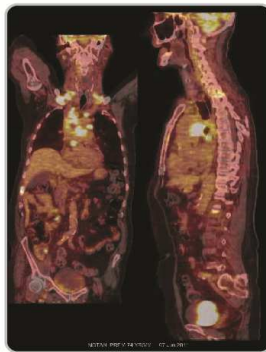

Cardiac CT



MRI

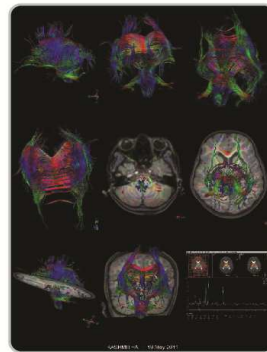

X-Ray



COLOR PRINTER
DS80


The essentials of imaging

Dye Sublimation Printing Technology

Ink ribbon : A thin polyester film coated with ink, in a sequence of yellow, magenta, Cyan, and Overcoat.

Paper : An exclusive receiving paper.

The dye ink is transferred to the paper, according to the amount of heat applied from the thermal Head.

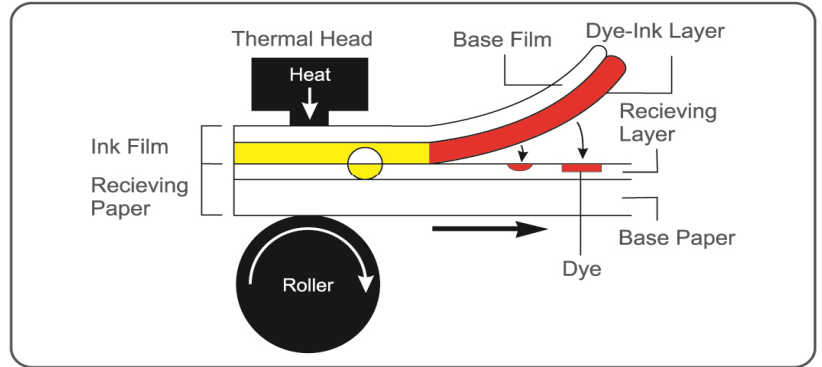


Dye Sublimation Printing Process


Various Applications for printing solutions DS 80 Medical Department That Use Printing Solutions

Specifications DS 80

Model	: Ds80
Print Method	: Dye Sublimation heat transfer method
Print Mode	: High-resolution mode
Resolution	: 600dpi
Print Size	: 8x12 203 x 305mm
Print Number	: 8x12 110sheets/ 1 Roll
Print Speed	: 8x12 46sec./sheet

Connectivity & Application

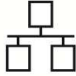LAN  DS 80 can be seamless connected to a lan network offering easy & fast connectivity to any imaging modality

USB  DS 80 can also be directly connected to any ultrasound or workstation PC using USB connection



KONICA MINOLTA

KONICA MINOLTA HEALTHCARE INDIA Pvt. Ltd.

Office No.515, 5th Floor, C-wing, 215-Atrium Centre Andheri (east), Mumbai 400 059. INDIA.
Tel. : + 91 - 22 - 61916969 Fax : + 91 - 22 - 61916996

URL: <http://konicaminolta.jp>